

syngo MammoReport

SP

System Manual

Disposal Instructions

© Siemens AG 2005

The reproduction, transmission or use of this document or its contents is not permitted without express written authority. Offenders will be liable for damages. All rights, including rights created by patent grant or registration of a utility model or design, are reserved.

Chapter	Page	Revision
All	All	01

Document revision level

The document corresponds to the version/revision level effective at the time of system delivery. Revisions to hardcopy documentation are not automatically distributed.

Please contact your local Siemens office to order current revision levels.

Disclaimer

The installation and service of equipment described herein is to be performed by qualified personnel who are employed by Siemens or one of its affiliates or who are otherwise authorized by Siemens or one of its affiliates to provide such services.

Assemblers and other persons who are not employed by or otherwise directly affiliated with or authorized by Siemens or one of its affiliates are directed to contact one of the local offices of Siemens or one of its affiliates before attempting installation or service procedures.

Copyright notice

The DICOM functionalities in the *syngo* MammoReport system are based on the libraries from Merge Technologies Inc.

Celsius is a trademark of Fujitsu Siemens Computers.

Eizo FlexScan L557 is a trademark of EIZO NANAO Corporation.

Windows XP Professional is a trademark of Microsoft Corporation.

SIMOMED SMM 21201P, DSB 2103-D-5MP, SCD 1897-M and *syngo* are trademarks of Siemens AG.

Coronis is a trademark of Barco.

Softcopy reading software by MeVis BreastCare,
copyright © 2001 - 2005 by MeVis BreastCare.



	Page
1 Overview of Hazardous Materials	1-1
General	1-1
Overview	1-1
Hazardous Material List	1-2
2 De-installation & Disposal	2-1
Disposal instructions	2-1
Hazardous materials.	2-1
Circuit board assembly, harnesses, monitors	2-1
Safety instructions.	2-1
Preparation	2-1
Hazardous Materials in the <i>syngo</i> MammoReport system	2-2
Arsenic / Arsenide and Cadmium	2-2
3 Changes to previous version	3-1

General

Regulations governing life-cycle management and refuse disposal hold the manufacturer/distributor responsible for the product. This includes environmentally responsible operation and disposal of the product.

The list of hazardous materials in these instructions gives an overview of the components requiring disposal. This information will assist the disposal department in locating these hazardous materials.

The disposal instruction provided in Chapter 2 should be followed when disposing of the system or replacing individual system parts.

Overview

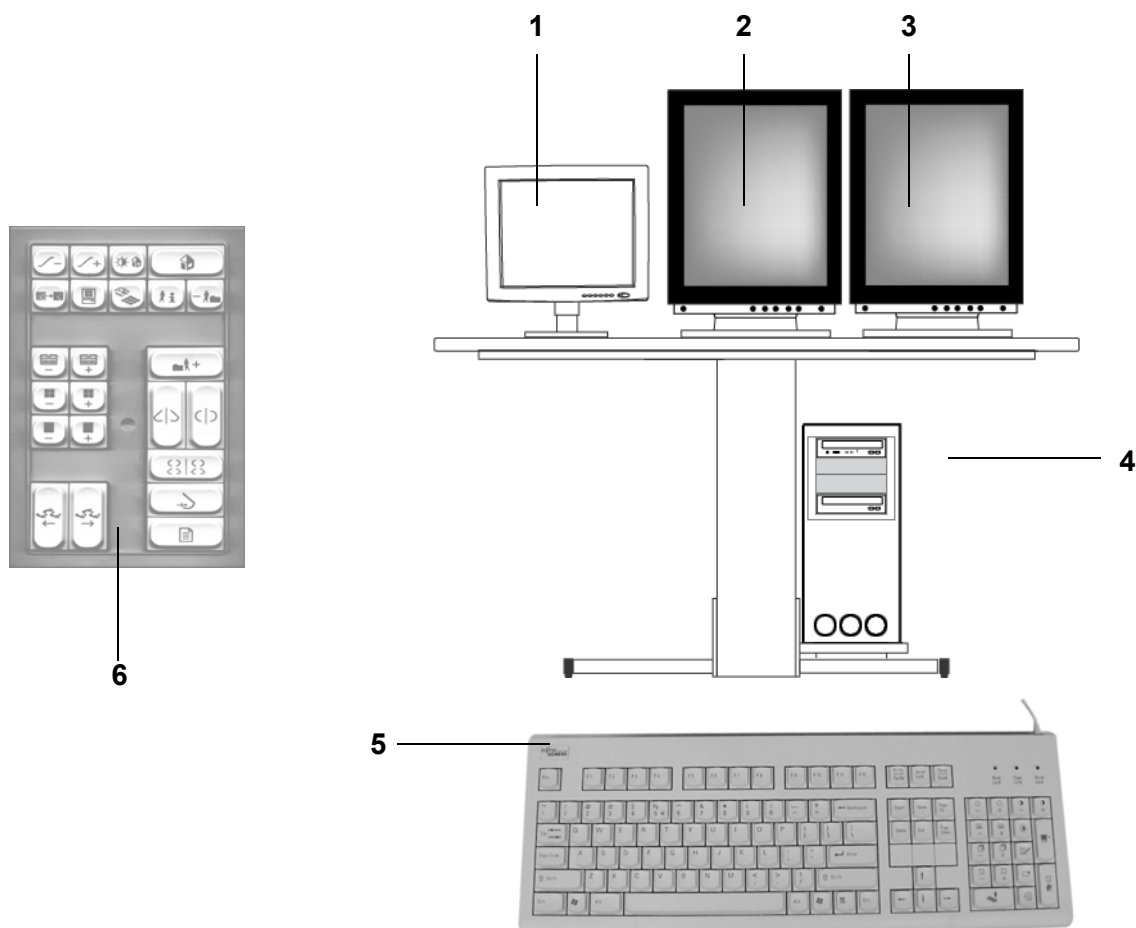


Fig. 1 Overview of syngo MammoReport

NOTE

The numbers in Fig. 1 are identical to the position numbers in the table on the next page.

Hazardous Material List

Pos.	Component	Hazardous Material	Quantity
1	19" LCD monitor/ syngo Monitor	arsenic / arsenide cadmium	<0,001% <0,001%
2	21" Monitor (CRT or TFT) / High Resolution Monitor	arsenic / arsenide cadmium	<0,001% <0,001%
3	21" Monitor (CRT or TFT) / High Resolution Monitor	arsenic / arsenide cadmium	<0,001% <0,001%
4	Workstation Celsius	arsenic / arsenide cadmium	<0,001% <0,001%
5	syngo Keyboard	arsenic / arsenide cadmium	<0,001% <0,001%
6	Workflow keypad	arsenic / arsenide cadmium	<0,001% <0,001%

Disposal instructions

Hazardous materials

Hazardous materials must be disposed of by an authorized waste disposal company according to industry standards. The hazardous materials appear in "bold/italic" in this chapter.

Circuit board assembly, harnesses, monitors

Circuit board assemblies, batteries, harnesses and monitors must be disposed of under supervision. Disposal procedures must be in accordance with industry standards.

Safety instructions

CAUTION

Note the general and product-specific safety instructions in the technical documents.

Preparation

WARNING

Risk of high voltage!

Disconnect the power cord of the workstation before starting the de-installation.

CAUTION

Remove the Lithium battery from the Computer.

Hazardous Materials in the *syngo* MammoReport system

Arsenic / Arsenide and Cadmium



Fig. 1 *syngo* MammoReport workstation

- Plastic parts of the workstation (PC) (1/Fig. 1):
 - Plastic parts containing **Arsenic / Arsenide and Cadmium**

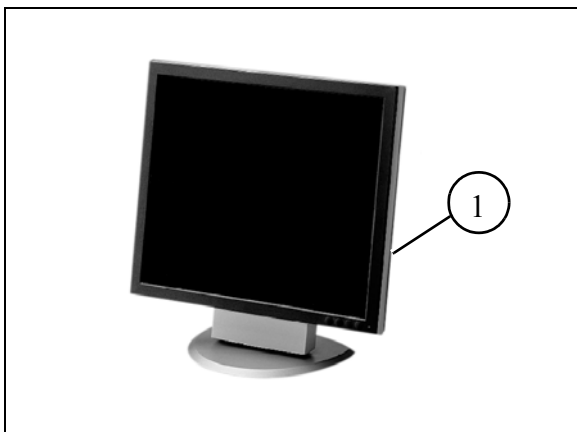


Fig. 2 19" syngo display

- Plastic parts of the 19" syngo display (1/Fig. 2):
 - Plastic parts containing **Arsenic / Arsenide and Cadmium**

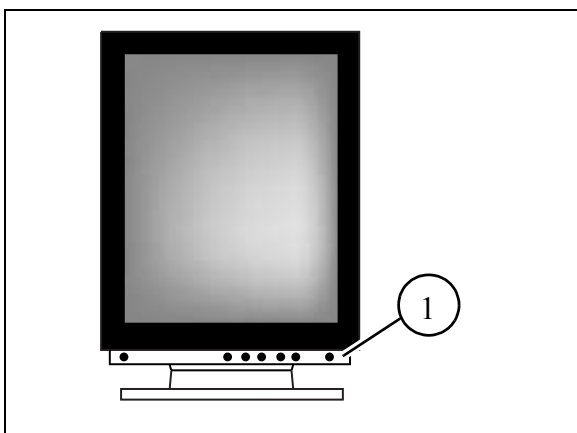


Fig. 3 21" High resolution monitor

- Plastic parts of the High resolution Monitor (CRT) (1/Fig. 3):
 - Plastic parts containing **Arsenic / Arsenide and Cadmium**

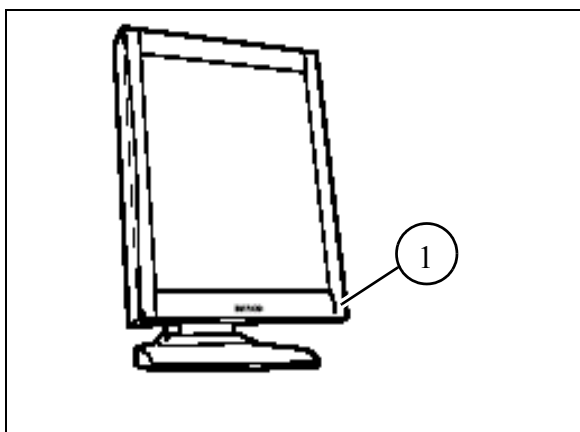


Fig. 4 21" Barco TFT monitor

- Plastic parts of the High resolution Monitor (Barco TFT)
(1/ Fig. 3):
 - Plastic parts containing **Arsenic / Arsenide and Cadmium**

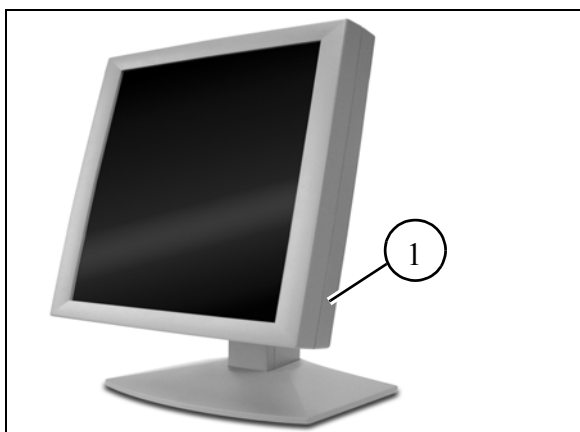


Fig. 5 21" Planar TFT monitor

- Plastic parts of the High resolution Monitor (Planar TFT)
(1/ Fig. 3):
 - Plastic parts containing **Arsenic / Arsenide and Cadmium**

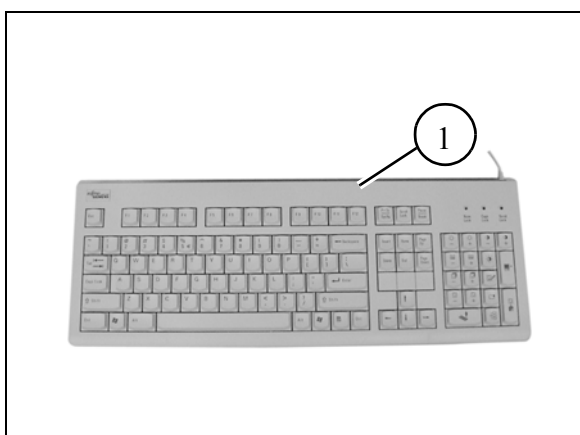


Fig. 6 syngo keyboard

- Plastic parts of the Keyboard
(1/ Fig. 6):
 - Plastic parts containing **Arsenic / Arsenide and Cadmium**

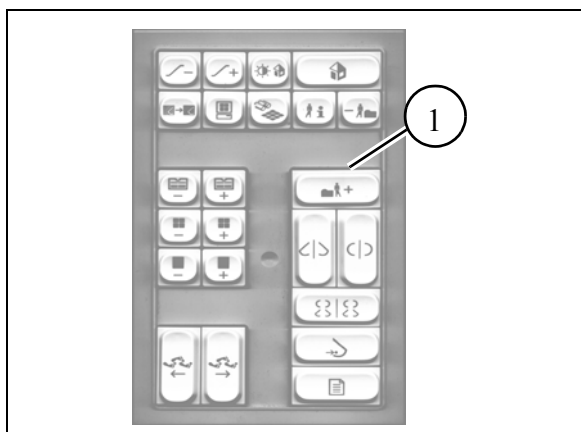


Fig. 7 Workflow keypad

- Plastic parts of the Workflow keypad (1/ Fig. 7):
 - Plastic parts containing **Arsenic / Arsenide and Cadmium**

Chapter	Page	Change
all		initial version

This page intentionally left blank.